

Technical Data Sheet


Alexa Fluor® 647 Rat anti-Histone H3 (pS28)**Product Information**

Material Number:	558217
Size:	50 tests
Vol. per Test:	20 µl
Clone:	HTA28
Immunogen:	Phosphorylated Human Histone H3 Peptide
Isotype:	Rat IgG2a, κ
Reactivity:	Tested: Human, Mouse, Reported: Cow, Drosophila, Hamster, Rat
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Histones are highly basic proteins that complex with DNA to form chromatin. Histone H3 is a ~15-kDa protein that is phosphorylated at serine 28 (S28), S10, and/or threonine 11 during mammalian cell mitosis and meiosis. The phosphorylation sites are located in the N-terminal tail, a region that is outside of the chromatin fiber and is thus accessible for interactions with agents that may regulate chromatin or specific gene activities. The phosphorylation states of the two serine sites during the cell cycle are highly regulated by Aurora B kinase and a PP1 phosphatase: S10 is in the phosphorylated state from late G2 phase to anaphase, while S28 is phosphorylated from prophase to anaphase. Furthermore, phosphorylation of histone H3 S28 may be mediated by other kinases in response to external stimuli. Evidence suggests that histone phosphorylation is involved in the regulation of chromosome condensation, cell division, and gene transcription.

The HTA28 monoclonal antibody reacts with histone H3 phosphorylated at S28 in its N-terminal tail. It does not recognize the unphosphorylated protein.

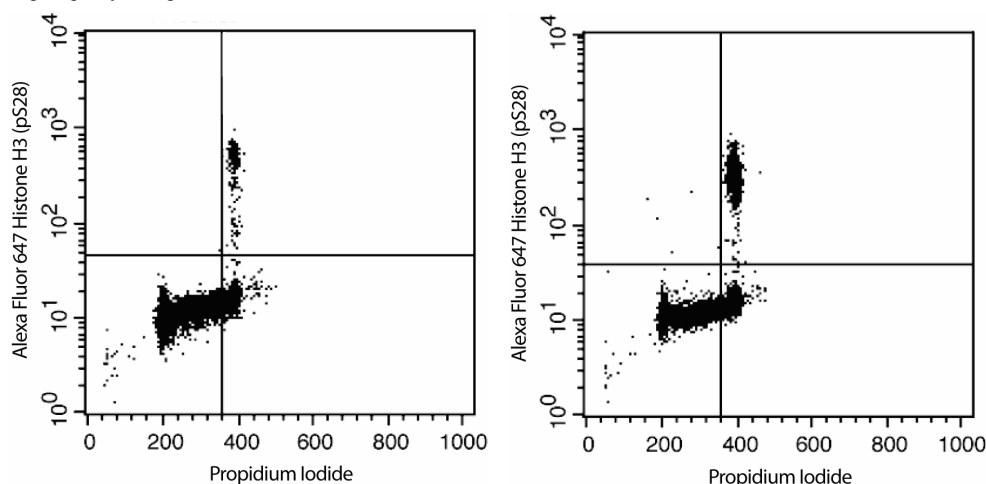

Analysis of Histone H3 (pS28) in proliferating mouse T lymphocytes. A log phase mouse cytotoxic T cell line was either untreated (left panel) or cultured with 1 µg/ml of the microtubule depolymerizer demecolcine (Sigma D7385) for 4 hrs (right panel). The cells were fixed and stained according to the Recommended Assay Procedure. Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system. In this example, the demecolcine treatment increased the frequency of G2/M phase cells from 20% to 52%, and the frequency of G2/M phase cells expressing Histone H3 (pS28) increased from 15% to 51% (upper right quadrants of the left and right panels, respectively).

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes**Application**

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

Recommended Assay Procedure:

Recommended Assay Procedure

1. Wash cell suspension twice with 1X PBS.
2. Fix the cells by adding ice-cold 70% ethanol drop-wise while vortexing the cell suspension, then storing them for at least 4 hours at -20°C in the 70% ethanol.
3. Aliquot ~1 million fixed cells per tube for staining. Wash them twice with 1X PBS, then once with stain buffer.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



4. Stain the cells with 20 μ l Alexa Fluor® 647 Rat anti-Histone H3 (pS28) in 80 μ l stain buffer for 20 minutes at room temperature, then wash them with stain buffer.
5. For optimum cell cycle analysis, the cells should be treated with RNase before staining with propidium iodide:
 - Treat the stained cells with 50 μ g RNase A (Sigma R5500) in 50 μ l 1X PBS for 30 minutes at 37°C. Without washing, stain DNA by adding 5 μ g Propidium Iodide (Sigma P4170) in 450 μ l staining buffer for at least 10 minutes at room temperature.
 Or
 - Stain the cells with 0.5 ml PI/RNase Staining Buffer (Cat. No. 550825) for 15 minutes at room temperature.
6. The cells are now ready for flow cytometric analysis.

Suggested Companion Products

<u>Catalog Number</u>	<u>Name</u>	<u>Size</u>	<u>Clone</u>
554656	Stain Buffer (FBS)	500 ml	(none)
550825	PI/RNase Staining Buffer	100 ml	(none)

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. Please refer to www.bdbiosciences.com/pharming/en/protocols for technical protocols.
3. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
4. This product is sold under license from Shigei Medical Research Institute, Okayama, Japan.
5. Alexa Fluor® 647 fluorochrome emission is collected at the same instrument settings as for allophycocyanin (APC).
6. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.
7. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
8. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
9. Alexa Fluor is a registered trademark of Molecular Probes, Inc., Eugene, OR.
10. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.

References

- Carmena M, Earnshaw WC. The cellular geography of aurora kinases. *Nat Rev Mol Cell Biol.* 2003; 4:842-854.(Biology)
- Cheung P, Allis CD, Sassone-Corsi P. Signaling to chromatin through histone modifications. *Cell.* 2000; 103:263-271.(Biology)
- Choi HS, Choi BY, Cho Y-Y, Zhu F, Bode AM, Dong Z. Phosphorylation of Ser28 in histone H3 mediated by mixed lineage kinase-like mitogen-activated protein triple kinase α . *J Biol Chem.* 2005; 280(14):13545-13553.(Biology)
- Dyson MH, Thomson S, Inagaki M, et al. MAP kinase-mediated phosphorylation of distinct pools of histone H3 at S10 or S28 via mitogen- and stress-activated kinase 1/2. *J Cell Sci.* 2005; 118(10):2247-2259.(Clone-specific)
- Furukawa K, Sugiyama S, Osouda S, et al. Barrier-to-autointegration factor plays crucial roles in cell cycle progression and nuclear organization in Drosophila. *J Cell Sci.* 2003; 116(18):3811-3823.(Clone-specific)
- Goto H, Tomono Y, Ajiro K, et al. Identification of a novel phosphorylation site on histone H3 coupled with mitotic chromosome condensation. *J Biol Chem.* 1999; 274(36):25543-25549.(Immunogen)
- Hirata A, Inada K, Tsukamoto T, et al. Characterization of a monoclonal antibody, HTA28, recognizing a histone H3 phosphorylation site as a useful marker of M-phase cells. *J Histochem Cytochem.* 2004; 52(11):1503-1509.(Clone-specific)
- Nowak SJ, Corces VG. Phosphorylation of histone H3: a balancing act between chromosome condensation and transcriptional activation. *Trends Genet.* 2004; 20(4):214-220.(Biology)
- Pascreau G, Arlot-Bonnemains Y, Prigent C. Phosphorylation of histone and histone-like proteins by aurora kinases during mitosis. *Prog Cell Cycle Res.* 2003; 5:369-374.(Biology)